

April 2002

PHOENIX MOUNT RO SYSTEM




We are pleased to announce the availability of a new mounting kit for attaching our Model PA2000 water purification system to the new Gambro Phoenix model dialysis machine.



Applied Water Engineering, Inc.
2759 South 300 West, Suite G
Salt Lake City, Utah 84115
Phone 801-486-7774, Fax 801-484-9265
Awater@xmission.com

The stainless steel mount bolts to the PA2000 and fits onto the Phoenix in place of the rear shelf. This transfers all the weight securely to the inner frame on the Phoenix. It is easy to remove and replace without tools. A small cinch strap holds the RO to the sturdy handle on the back of the Phoenix.

The kit includes bottle mounts and quart bottles for Actril™ and bleach. Another bracket has holes for relocating the disinfection wands that come with the Phoenix.

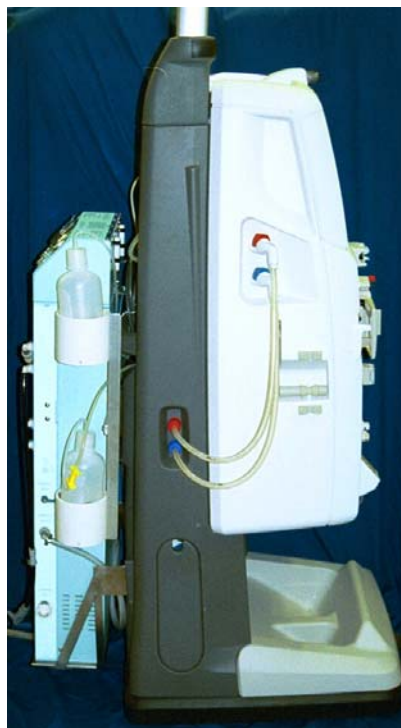

The PA2000 is now in the fifth year of production. It is unique among RO systems in its high rated output, differential pressure gauges for both prefilters, triple pressure regulators, and full fit RO module with Circsert™.



Applied Water Engineering, Inc.
2759 South 300 West, Suite G
Salt Lake City, Utah 84115
Phone 801-486-7774, Fax 801-484-9265
Awewater@xmission.com